

MS CEILING MS CEILING PLUS

Minor surgery and advanced exploration luminaire designed to be ceiling mounted. The lamp incorporates a LED light source and a highly precise articulated arm that allows the luminaire to be adjusted in any position. The quality of the materials used in its manufacture makes it especially suitable for intensive use in surgery and examinations.

Exceptional quality of LED light

The powerful LED light source offers exceptional light quality and a lifespan of more than 50,000 hours.

The high performance of the colour Ra 97% depicts the exceptional quality of white light for an excellent distinction of fabrics and its reproducing all visible colours, naturally and precisely.

Ergonomic design

The circular shape of the head and its handle in the lower part facilitate its exact positioning, fulfilling the ergonomic requirements for any use.

Designed with clean and soft shapes which facilitate hygiene and cleanliness in medical environments.

Ceiling Support

The ceiling support is included. It is manufactured specifically for the MS CEILING family lamp.

The MS CEILING family lamp is not compatible with any other accessory.

It is designed to be installed in ceilings up to 3 meters.

Versions

MS CEILING	45.000 lux at 50 cm On/Off version
MS CEILING PLUS	65.000 lux at 50 cm Adjustable light intensity version



MS CEILING PLUS

Technical Details	(REF. M187001C)
Light source	LED 12W
Lux at 50 centimeters	65.000 lux
Adjustable light intensity	Yes
Colour temperature	4.500 °K
Colour rendering	> 97%
Lamp lifetime	50.000 hours
Arm length	105 cm
Ceiling support length	100 cm maximum
Colour	White

Medguard Professional Healthcare


Tel: 01 835 2378

Email: orders@medguard.ie

Fax: 01 969 5009

Web: www.medguard.ie

Unit 6/7, Block 4, Ashbourne Business Park, Ashbourne, Co. Meath, Ireland

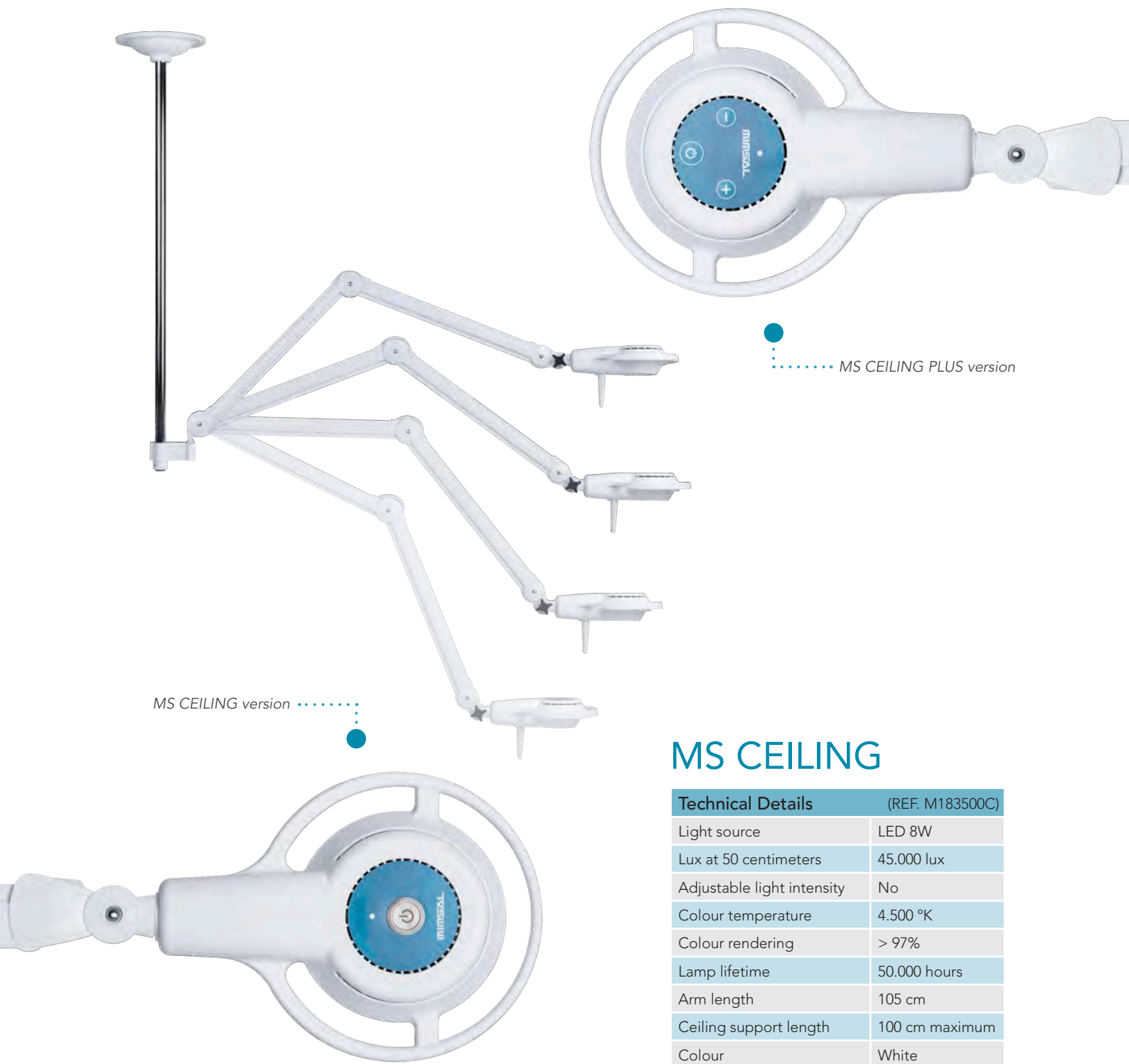

MS CEILING version

MS CEILING PLUS version

MS CEILING

Technical Details	(REF. M183500C)
Light source	LED 8W
Lux at 50 centimeters	45.000 lux
Adjustable light intensity	No
Colour temperature	4.500 °K
Colour rendering	> 97%
Lamp lifetime	50.000 hours
Arm length	105 cm
Ceiling support length	100 cm maximum
Colour	White

