

Premier Detectors X-Ray Detectable Swabs



The 365 Healthcare range of Premier Detector swabs are carefully made to the highest levels of quality and safety.

Detector X-Ray detectable gauze swabs are manufactured to EN14079 standards and are recognised for their quality, absorbency and strength.

Cotton Gauze

- Soft and highly absorbent gauze to increase efficiency.
- Consistent quality of gauze construction, absorbency levels and weight.

X-Ray Detectability

- Higher than EU standard barium sulphate content (60-62%) ensuring maximum visibility on all commonly used x-ray devices.
- For maximum patient safety, the x-ray detectable strips are securely attached using heat and pressure providing consistent adherence throughout the body of the gauze.
- Metal detection testing employed during manufacturing process to ensure patient safety.

Taped Swabs

- A consistent stitch pattern with the tape end embedded 5 cm into the swab for a strong and secure attachment.
- Tape presented on the outer surface of the swab for speedy counting and minimal swab handling.


Presentation

- Length and position of the red tie cord ensures easy removal from swab bundle.
- Quality controlled at three stages during the manufacturing and packing processes to ensure correct bundles of 5.
- Double wrapped in line with aseptic opening practices to minimise contamination risks.

Peel Open Pouches



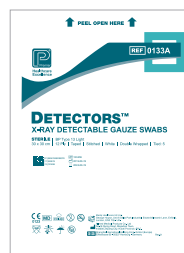

- Large opening lip enables easy opening and safe presentation of sterile product.practice to minimise contamination risks.

- Manufactured to BP/EN 14079 Standard
- Extensive choice of sizes and ply's to satisfy most requirements from one supplier
- Extensive non sterile range also available



Materials:	100% Cotton Gauze Type 13
Colours	White and Green
Intended Use:	For intra operative use only
Inner Wrap:	Plain white paper (sterile, double wrapped only)
Outer Wrap:	Thermoform peel open pouch
Sterilisation:	ETO

PREMIER		TIED 5'S DOUBLE WRAPPED   STERILE 							
CODE	PLY	WIDTH	LENGTH	COLOUR	STERILE	STITCHED	TAPED	INNER QTY (TIED 5'S)	OUTER QTY (TIED 5'S)
0102	8 ply	7.5 cm	7.5 cm	Green	●			80	160
0103	12 ply	7.5 cm	7.5 cm	White	●			80	160
0150	12 ply	10 cm	10 cm	Green	●			70	140
0106	12 ply	10 cm	10 cm	White	●			70	140
0130A	12 ply	22.5 cm	22.5 cm	White	●	●	●	10	80
0160	12 ply	22.5 cm	22.5 cm	White	●	●		24	48
0133A	12 ply	30 cm	30 cm	White	●	●	●	10	40
0166A	12 ply	30 cm	30 cm	White	●	●		10	40
0165	12 ply	30 cm	30 cm	White	●			14	28
0142A	12 ply	45 cm	45 cm	White	●	●	●	8	24
0170	12 ply	45 cm	45 cm	White	●	●		6	12
0107A	16 ply	10 cm	10 cm	White	●			80	160
0113A	16 ply	15 cm	2.5 cm	White	●	●		80	400
0122A	16 ply	36 cm	11 cm	White	●			15	60
0116A	24 ply	15 cm	10 cm	White	●	●		15	90
0136A	24 ply	30 cm	30 cm	White	●			6	24
0101	32 ply	5 cm	5 cm	White	●	●	●	80	160
0104A	32 ply	10 cm	7.5 cm	White	●			15	120
0110A	32 ply	10 cm	10 cm	White	●			15	120

PREMIER		TIED 5'S BULK WRAPPED    							
CODE	PLY	WIDTH	LENGTH	COLOUR	NON-STERILE	STITCHED	TAPED	INNER QTY (TIED 5'S)	OUTER QTY (TIED 5'S)
7139	6 ply	45 cm	45 cm	White	●		●	10	120
7103	8 ply	7.5 cm	7.5 cm	White	●			20	1,000
7012A	8 ply	7.5 cm	7.5 cm	Green	●			100	6,000
7109	12 ply	10 cm	10 cm	White	●			20	1,000
7156A	12 ply	10 cm	10 cm	Green	●			20	800
7129	12 ply	22.5 cm	15 cm	White	●			20	160
7130	12 ply	22.5 cm	22.5 cm	White	●	●	●	20	160
7131	12 ply	22.5 cm	22.5 cm	White	●			20	160
7133	12 ply	30 cm	30 cm	White	●	●	●	10	80
7132	12 ply	30 cm	30 cm	White	●			10	80
7143	12 ply	45 cm	45 cm	White	●			10	40
7104	12 ply	7.5 cm	7.5 cm	White	●			20	1,000
7110	16 ply	10 cm	10 cm	White	●			20	480
7125	16 ply	36 cm	11 cm	White	●			20	160
7116	16 ply	15 cm	2.5 cm	White	●	●		20	240
7122	16 ply	25 cm	2.5 cm	White	●	●		20	240
7134	18 ply	30 cm	30 cm	White	●		●	10	40
7119	24 ply	15 cm	10 cm	White	●			20	240
7136	24 ply	30 cm	30 cm	White	●		●	10	40
7101	32 ply	5 cm	5 cm	White	●			20	1,000
7107A	32 ply	10 cm	7.5 cm	White	●			10	400
7013A	32 ply	10 cm	7.5 cm	Green	●			50	2,000
7113	32 ply	10 cm	10 cm	White	●			20	240