


BD VENFLON PRO SAFETY IV CATHETERS 20G 32 mm - 393224 - sterile

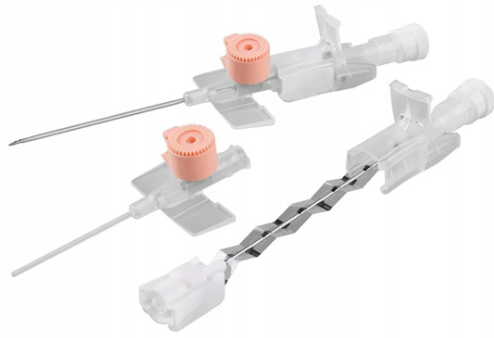

Product Code	23713
Unit of sale	box of 50 pcs
Minimum order	1
Type	Medical device
Class	II A

RDM (NSIS)	1883246
CND	C0101010201
GMDN	64574

Description

BD VENFLON PRO SAFETY NEEDLES - sterile

- BD code: 393224
- Gauge: 20g
- Gauge mm: 1,1x32 mm
- Colour: pink

BD Venflon™ Pro Safety Needle Protected IV Cannula is a ported, passive safety peripheral IV catheter (PIVC) designed to help minimise the risk of unanticipated blood spatter and needlestick injuries (NSIs). BD Vialon™ catheter material (polyurethane, radiopaque and biocompatible) softens in the vein up to 70%, enabling longer dwell times, minimizing restarts and reducing the risk of mechanical phlebitis. Ergonomic design to accommodate a variety of techniques. Back-cut needle design to offer optimal skin and vein penetration. Flexible wings help provide comfortable securement for the patient. "Snap" cap helps prevent port being accidentally opened. Removable white vent cap provides flexibility to address a variety of needs.

Multilingual instructions in: GB, FR, IT, ES, PT, DE, NL, PL, DK, FI, SE, NO, EE, LV, LT, CZ, SI, SK, GR, BG, HR, HU, RO, UA, TR, RU, KZ, Arabic.